

Base 2™ Outrigger Components

Intended Purpose

These non-sterile components are used to create a custom, lightweight outrigger support for static or dynamic splinting of hand and wrist.

Indications

Optimal positioning of the wrist and digits post injury or diagnosis using a combination of the outrigger components.

Contraindications


Skin irritation. Allergy to Nylon, Aluminum, Suede.

Instructions for Use / Proper Fit

All components are interchangeable and can be used to customize multiple outrigger splints.


Note: Continuously monitor the patient's condition and the splinting site during use. Make adjustments as necessary. Do not use excessive pressure on bony prominences or sensitive areas during application. Advise patient to monitor for adverse reactions or complications that may arise during use, such as swelling, numbness, or changes in skin color.

All aluminum components may be custom contoured by bending with tools or by hand, over a table edge, or using a table vise.

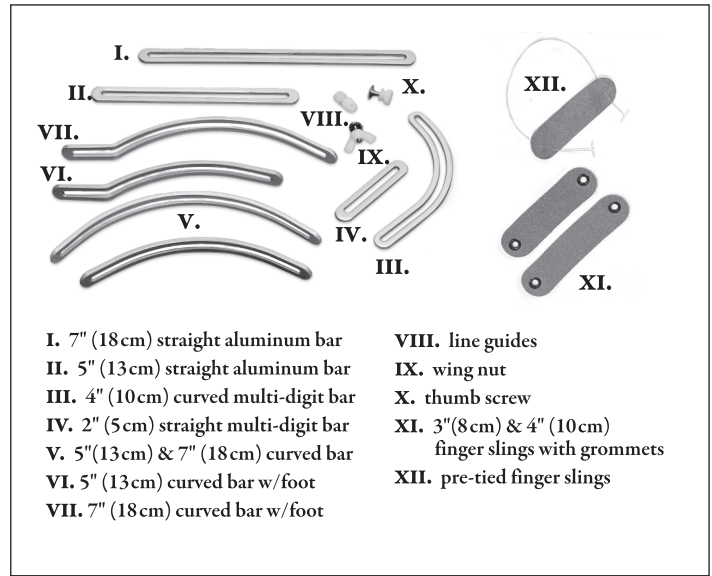

The following photos show examples of how the components may be used in combination to attach an outrigger to a thermoplastic splint.


Single Finger Extension Splint

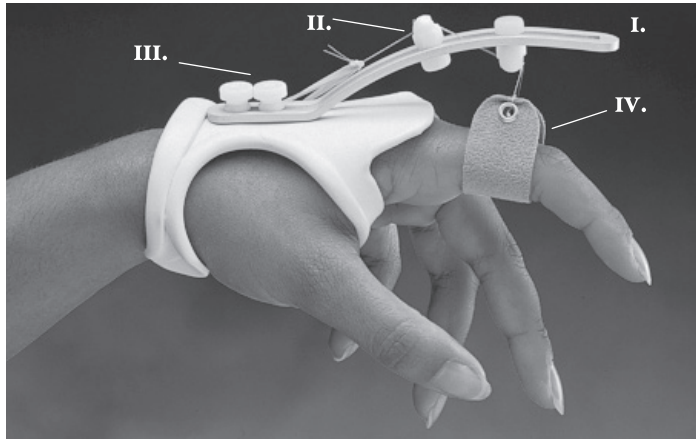

Parts Used:

- I. 5" (13cm) curved bar with foot
- II. line guides
- III. thumb screws
- IV. finger sling with grommets

Multi-Digit MP Extension Splint

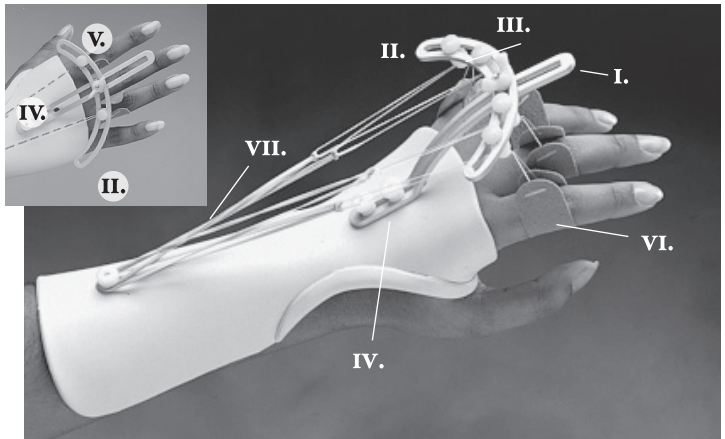

Parts Used:

- I. 5" (13cm) curved bar with foot
- II. 4" (10cm) curved multi-digit bar
- III. line guides
- IV. thumb screws
- V. wing nut
- VI. pre-tied finger slings
- VII. rubber bands

Instructions For Care: Hand wash in cool water and mild soap. Air dry.

Warning:

- Avoid direct contact of components with open wounds or damaged skin.
- Do not attempt to punch or drill holes in the thermoplastic while it is on the patient's hand.
- Choking Hazard - Keep out of reach of children.
- Discontinue use or adjust fit if patient shows signs of irritation, impaired circulation, increased pain or discomfort, or allergic reactions such as redness, itching, tingling, rash, or color changes in the affected area.
- Use of the Outrigger Kit for other than its intended purposes may cause injury.

Caution:

- Patient should have appropriate skin hygiene practices to prevent complications associated with prolonged splint use.
- Store Outrigger Kit and components in a clean and dry area.

60°F (16°C) 85°F (29°C)

UK CA CE

EC REP UK REP

CHN/156

MD

CH REP

North Coast Medical, Inc.
135 E Main Ave, Suite 110
Morgan Hill, CA 95037-U.S.A.

ICON (LR) Limited
South County Business Park
Leopardstown, Dublin 18
D18 X5R3, Ireland

Allenspach Medical AG
Dornacherstrasse 5
4710 Balsthal, Switzerland



North Coast Medical, Inc.
www.ncmedical.com