



NEW *products!*

VOL. 1.1

PPE Supplies

- Face Mask
- Safety Glasses & Face Shield
- Hand Sanitizers & Surface Disinfectant

Orthotics

- Progress™ Ultra-Ortho
- Comfort Cool® Gladiator™
Wrist Orthotic
- Ortho Foam™ High-Force
PIP Extension Orthotic
- Vasyli™ Foot Orthotics
- Norco® Nite-Nite™ Wrist Support

Therapy Products

- Stretching
- Balance & Self-Massage
- Cold Packs
- E-Stim Electrodes
- Dry Needles



FREE Norco® Stretch Strap!


Add a Norco® Stretch Strap (NC84527) to any order. Use promo code **NP21S** at checkout. Good through 05/31/21

MEDICAL GRADE

DISPOSABLE FACE MASKS



- **ASTM F2100 Level 2**
- 3-ply design provides BFE and PFE > 95%.
- Adjustable **nose clip** significantly reduces air penetration.
- FDA registered manufacturer.

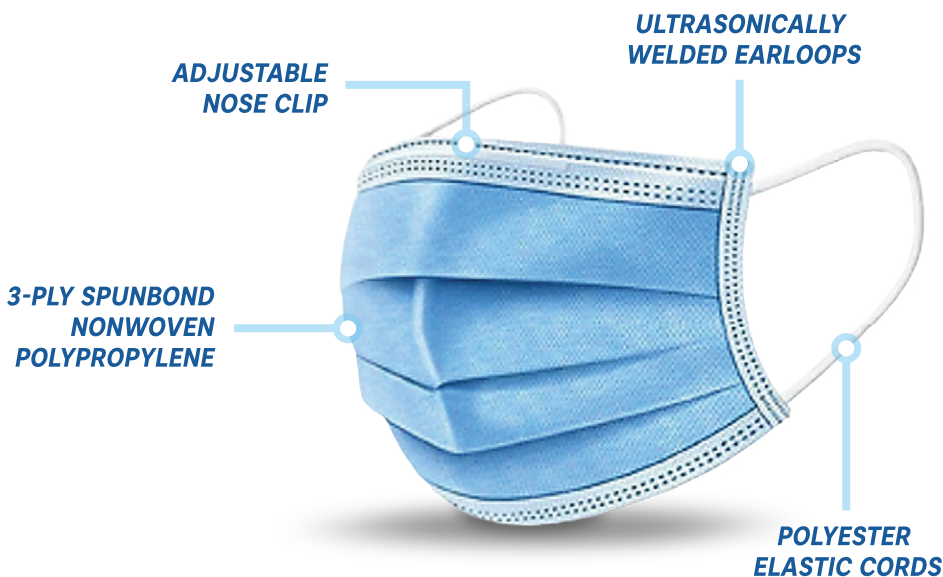

Dimensions: 7 x 3.75" (17.7 x 9.5 cm)
Not made with natural rubber latex.

NC36651 Box of 50

**EVERYDAY
EVERYONE™
PRICING**

AFFORDABLE & DISPOSABLE

To order, log in at **www.ncmedical.com**

SAFETY GLASSES FOR ANY WORK PLACE

OTG Safety Glasses

- Ideal for applications that require eye protection.
- H2X anti-fog, water and scratch-resistant coating.
- Fits over prescription glasses.
- Made of dielectric material.
- Polycarbonate lens.

*Temples extend
from browbar to
avoid interference
with prescription
frame temples.*



Lens Specifications

Visible light transmitted 380-780 nm = 96%

UVA, UVB blocked 200-380 nm = 100%

Blue light blocked 400-700 nm = 9%

IR blocked 780-2000 nm = 18%

NC28846 Box of 12

To order, log in at www.ncmedical.com



SINGLE-USE FACE SHIELD

Protection for the entire face.

Lightweight, disposable, cost effective product that provides comfortable coverage.

- Made with optically clear, 7 mil thick film.
- Full length 9" (23cm) shield protects face from direct splatter.
- Elastic webbing for comfortable, yet secure fit.
- Not made with natural rubber latex.
- Made in USA.



SINGLE-USE FACE SHIELD

Pack of 10
NC36652

To order, log in at
www.ncmedical.com

WAGING WAR ON GERMS & BACTERIA



GERM WAR

HAND SANITIZER

**FIRST LINE OF DEFENSE
AGAINST COMMON COLD AND FLU**

Active Ingredient: Ethanol 62%



HELPS FIGHT BACTERIA AND VIRUSES

**WHEN USED PROPERLY,
HELPS CLEAN AND SANITIZE HANDS**

**PACKAGED FOR EASE OF USE
AND TRAVEL MOBILITY**

MOISTURIZING WITH ALOE

**NC17379-1
16.9 oz Pump**

**EVERYDAY
EVERYONE™
PRICING**

To order, log in at **www.ncmedical.com**

Safetec
SaniZide^{PRO}
Ready-to-Use Liquid Disinfectant

Effectively KILLS COVID-19,
 Influenza, Norovirus, MRSA,
 Poliovirus *and more*

1-STEP 1-MINUTE

Hospital-grade, alcohol-based formula
 kills 49 microorganisms in 1 minute

- Ideal for use in hospitals, nursing homes, medical or dental offices, ambulances, pharmacies and other healthcare related institutions.
- EPA registered surface disinfectant that combines cleaning and disinfecting in one simple step with its ready to use, alcohol based hospital grade formula.
- Safe for repeated use on hard, non-porous surfaces.

EPA Registration number 88494-3

NC34821 Spray Bottle
 32 oz. (907 g)
 6 bottles

**EVERYDAY
 EVERYONE™
 PRICING**



To order, log in at **www.ncmedical.com**

HAND SANITIZING SOLUTIONS

for a healthier you!

giovanni®
SKIN & BODY CARE

Moisturize and reduce germs and bacteria instantly with effective and high-performance hand sanitizing solutions.


Powerful, fragrance-free Hand Sanitizer is infused with nourishing antioxidant Vitamin E, Chamomile, and certified organic Aloe Vera.

Mixed Hand Sanitizing Towlettes send germs packing while soothing hands with certified organic oils of Grapefruit, Lavender, Peppermint and Spearmint. Available in conveniently packed travel sizes.



Refresh With Aloe Vera

A gentle solution for non-greasy skin care — hydrates and soothes while providing powerful antioxidants to help rejuvenate each pore.

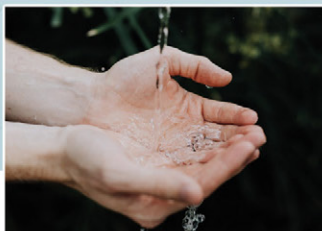

Fight Bacteria With Ethyl Alcohol

The powerful 70% Ethyl Alcohol in our hand sanitizer gel helps eliminate germs — perfect for cleansing and refreshing at any time.



Nourish With Botanical Extracts

These cruelty-free hand sanitizers are infused with nutrient-rich botanical extracts to promote cleaner, healthier, and more vibrant skin.

Clean and Refresh Skin, Anytime, Anywhere!

Giovanni® Antibacterial Hand Sanitizer, Fragrance-free Volume Quantity

NC17386	8.5 oz (250 ml) Pump	1
NC17387	24 oz (710 ml) Pump	1
NC17385-24	2 oz (59 ml) Bottle	24

Giovanni® Antibacterial Towlettes, Scented Scent Quantity

NC17388	Mixed	1 pack (20 individually wrapped towlettes)
---------	-------	--

To order, log in at www.ncmedical.com

ULTRA- ORTHO™

Multipurpose Orthotic



Provides wrist positioning and adjustability for both flexion and extension of the hand, wrist, thumb and digits.

Allows Metacarpophalangeal (MCP) joints to be placed into **extension** for extensor tendon injuries, or **flexion** for joint injuries.



C-BAR
Maintains webspace of thumb allowing a good sweeping position for adjustable abduction, extension and flexion.



"The Progress Ultra-Ortho is a functional solution for therapists and physicians. The orthotic is versatile and customizable for a wide variety of diagnoses, both orthopedic and neurological. The skilled craftsmanship and superior materials provide confidence knowing patients are getting a product that will therapeutically protect, progress, and meet all functional goals."

- Greg Pitts, OTD, OTR/L, CHT



The innovative design combines a malleable frame with dense foam padding for a lightweight, soft orthosis that supports like a moldable plastic splint.

**ULTRA-
ORTHO™**
Multipurpose Orthotic



Independent thumb strap allows custom adjustability for securing and supporting thumb in place.

All straps are fastened with a rivet to meet L-Code requirements. Each strap is rounded to prevent shearing and offer maximum comfort and control.



Ultra-Ortho™ is biomechanically designed to immobilize and offer adjustable support to meet a variety of therapeutic goals.

The orthotic is easily adjustable and does not cause any impingement or pressure areas when straps are fastened. This allows the orthotic to be used immediately for a variety of different types of pathologies, maximizing the functional outcome and progress as needed without refitting for a new orthotic.

Applications include but are not limited to: contractures due to Dupuytren, head injury, stroke, spinal cord injuries, CVA, arthritic conditions or burns.

Direction of strapping biomechanically controls the hand, wrist, and forearm.

Left	Right	Size
NC25560	NC25561	X-Small
NC25562	NC25563	Small
NC25564	NC25565	Medium
NC25566	NC25567	Large

One model, multiple applications

The versatile Ultra-Ortho™ eliminates the need of numerous other orthoses, including these Progress™ orthotics:



Progress™
I.F. (Individual Finger)
Ball Orthosis



Progress™
Neutral Resting
Orthosis



Progress™
Neutral Hand Cone
Orthosis



Progress™
Dorsal Anti-Spasticity
Orthosis



Progress™
Palmar Bar Wrist
Orthosis

To order, log in at www.ncmedical.com



GLADIATOR™

Superior Grade Materials

Premium Italian Fabrics
for style and comfort



Anti-Slip Beads
for traction and a secure fit



Wrist Orthosis

Biomechanical Advantage

- Counterforce straps
- Formfitting
- 2/3 forearm length
- Removable stays
- Breathable
- Lightweight

Perforated Ariaprene®
provides comfortable compression and breathability



D-Rings
anodized aluminum for added strength

Low Profile Hook
injection molded featuring easy pull tabs



"Gladiator is designed to comfortably protect post-surgical and chronic wrist disorders. The unique bio-mechanical design allows the physician to provide a custom fit while maximizing occupational function."

- Greg Pitts, OTD, OTR/L, CHT

Comfort Cool® Gladiator™ Wrist Orthosis

	Size	Wrist Circumference	Hand
NC86660	Small	6" to 7" (15 to 18 cm)	Left
NC86661	Small	6" to 7" (15 to 18 cm)	Right
NC86662	Medium	7" to 8" (18 to 20 cm)	Left
NC86663	Medium	7" to 8" (18 to 20 cm)	Right
NC86664	Large	8" to 9" (20 to 23 cm)	Left
NC86665	Large	8" to 9" (20 to 23 cm)	Right
NC86666	X-Large	9" to 10" (23 to 25 cm)	Left
NC86667	X-Large	9" to 10" (23 to 25 cm)	Right

Comfort Cool® Gladiator™ Wrist & Thumb Orthosis

	Size	Wrist Circumference	Hand
NC86650	Small	6" to 7" (15 to 18 cm)	Left
NC86651	Small	6" to 7" (15 to 18 cm)	Right
NC86652	Medium	7" to 8" (18 to 20 cm)	Left
NC86653	Medium	7" to 8" (18 to 20 cm)	Right
NC86654	Large	8" to 9" (20 to 23 cm)	Left
NC86655	Large	8" to 9" (20 to 23 cm)	Right
NC86656	X-Large	9" to 10" (23 to 25 cm)	Left
NC86657	X-Large	9" to 10" (23 to 25 cm)	Right



Distal Strap

extra-long, customizable strap can be cut to specific length for optimal support. Two-sided hook fastener attaches easily to strap.

Bicep Cutout
allows for active range of motion



To order, log in at
www.ncmedical.com

ORTHO FOAM™

High Quality Economical Price

FINGER ORTHOTICS



High-Force

PIP Extension Orthotic

PIP joint and tendon correction.



- Use to correct PIP extension lag, prevent volar plate stiffness, and protect injured extensor tendon.
- Contraindicated for fixed contractures, patients with fragile skin or severe edema.
- Heavy duty (0.047) black, stainless-steel wire construction. Easy to fit and bend wires to adjust extension force.
- Large closed-cell, contoured foam pads for improved fit.
- Not made with natural rubber latex.



"The Ortho Foam High-Force PIP Extension Orthotic is the only product I've found that provides a superior biomechanics design that protects healing tissue, and maximizes patient comfort and functional outcome."

– Greg Pitts, OTD, OTR/L, CHT

Ortho Foam™ High-Force-PIP Extension Orthotic

	Size	DPC to DIP Length
NC85423-1	X-Small (AA)	2¼" (5.7 cm)
NC85423-2	Small (A)	2⅝" (6.7 cm)
NC85423-3	Medium (B)	3" (7.6 cm)
NC85423-4	Large (C)	3½" (8.9 cm)
NC85423-5	X-Large (D)	4" (10 cm)

To size, measure from the distal palmar crease to the distal finger crease (finger extended).




To order, log in at **www.ncmedical.com**

Foam Arm Support with Vinyl Coating

**Use to help reduce edema after surgery
or protect flaccid upper extremity
after stroke.**

- Contoured, dense arm support for positioning the injured or weak hand comfortably and safely near the body, or elevating the swollen hand above the heart to aid in fluid return.
- Holds the hand at a comfortable elevation during treatment and exercise.
- Fits on a lapboard, table or bed without slipping.
- Maintains the hand 5" (13 cm) above tabletop.
- Measures 20" in length and 5 1/2" wide (51 x 14 cm).
- Made of solid, open-cell foam.
- Features a contoured trough that supports the forearm.

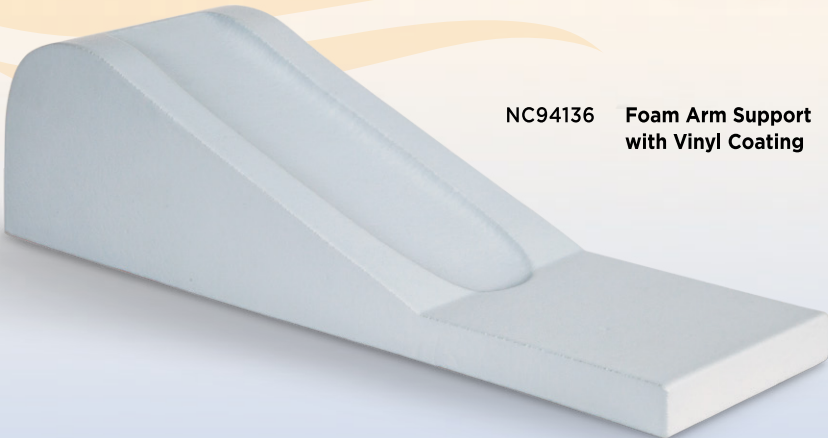

NC94136 **Foam Arm Support
with Vinyl Coating**




Smart Splint Pan




Sterilizes **99.99%**
of microbes using
high heat and UV light

- Touch screen user interface with parameter adjustability, thermoplastic sheet presets, user presets, auto on/off, education library and more.
- Transparent lid with accordion fold for visibility, heat maintenance and easy access to thermoplastic sheets.
- Heat cold water to precise thermoplastic temperature in less than 20 minutes.
- Automatic drainage pump for quick and easy water evacuation. A built-in water level sensor turns off the pan when water level is too low.
- Thermoset composite construction for durability, insulation, flammability resistance and aesthetic appeal.
- Recessed thermostat for safe and even heating.
- External Dimensions: 20" x 23" x 9" (51 x 58 x 23 cm).
- Tub Dimensions: 13" x 19" x 6" (33 x 49 x 15 cm).
- Glass height (fully raised) measures 9.5" (24 cm). Total height reaches 18.5" (47 cm).
- IEC61010-2 safety rating for hospital and international compliance.

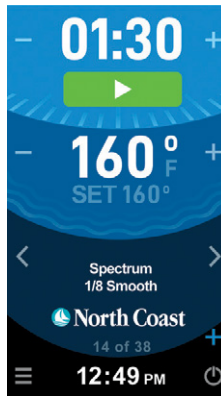

Home screen with factory presets of thermoplastic sheet recommendations.

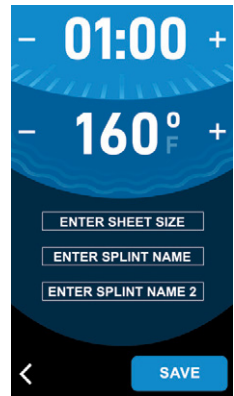

Create custom presets to fit clinician preferences and non-programmed thermoplastics.

**Strava Smart Splint Pan
NC14599**

MADE IN THE USA
BY **VETERANS**

To order, log in at www.ncmedical.com

Preferred® Black



Our new, most versatile, non-stick splinting material can be used by novice or experienced users alike.

Resistance to Stretch

- Minimum to moderate resistance to stretch.
- Tolerates minimum to moderate handling during splint molding.
- Handle on a horizontal plane and form splint in a gravity-assisted position.
- Gentle, smooth motions yield superior results.

Drape

- **Moderate drape** with only light to moderate pressure.
- Contours are quick to form.

Memory

- **Minimum to moderate memory.**
- Splints can be modified a few times without thinning.
- Requires less handling and maintains contours while cooling.

Rigidity

- Molded splints have a **moderate rigidity**.
- A circumferential splint is flexible enough for easy donning and doffing, yet strong enough for immobilization.

Bonding

- The **NS coating** gives a slightly tacky feel that helps hold the material in place, without the material pulling or falling away.
- A wet bond will hold but can be pulled apart once cooled.
- Dry heating two surfaces will form a strong bond. For a permanent bond, remove coating with an emery cloth and apply dry heat, then press together firmly. Solvent not required.

Surface Finish

- **Moderate resistance to fingerprints** and resists showing stains.
- **Smooth, self-sealing** cut edges remain sealed even after reheating.

Applications

- **1/8" Preferred® Black** can be used for almost any splinting application. Use for dynamic and static splinting; for adult and pediatric positioning; and for post-burn and post-surgical care.
- **3/32" Preferred® Black** is thin and lightweight, yet strong, and can be used for forearm and hand-based splints. Ideal for arthritis splints, post-operative splints, and for pediatric splints.

Working Time (1/8")

- Allows three to five minutes of molding time after softening in 160°F (71°C) water.

Not made with natural rubber latex.

1/8" (3.2 mm) Preferred® Black Smooth

NC12813	18" x 24" (46 x 61 cm)	Case (4)
NC12813-1	18" x 24" (46 x 61 cm)	Sheet (1)

1/8" (3.2 mm) Preferred® Black Perforated 1%

NC12814	18" x 24" (46 x 61 cm)	Case (4)
NC12814-1	18" x 24" (46 x 61 cm)	Sheet (1)

3/32" (2.4 mm) Preferred® Black Smooth

NC12815	18" x 24" (46 x 61 cm)	Case (4)
NC12815-1	18" x 24" (46 x 61 cm)	Sheet (1)

3/32" (2.4 mm) Preferred® Black Perforated 1%


NC12816	18" x 24" (46 x 61 cm)	Case (4)
NC12816-1	18" x 24" (46 x 61 cm)	Sheet (1)

Thickness Guide

1/8"	3.2 mm
3/32"	2.4 mm

Perforation Guide

- Smooth
- Perforated, 1%

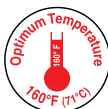

Performance characteristics similar to TailorSplint® Black

TailorSplint is a registered trademark of Performance Health Holdings, Inc.

To order, log in at **www.ncmedical.com**



NEW COLORS GREAT QUALITY SAFETY-SEALED CONTAINERS



Rainbow™ EXERCISE PUTTY

**Promote hand-muscle strength
and increase endurance.**



- Exercising with Rainbow™ Exercise Putty helps promote increased range of motion, finger-flexion and tendon-gliding.
- Rainbow™ color-coded system associates resistance levels with color.
- The smaller portion cups are perfect for most small and medium sized hands.
- Composed of silicon rubber.
- Not made with natural rubber latex.



Soft Orange

NC98200-02	2 oz. (57 g)
NC98200-03	3 oz. (85 g)
NC98200-04	4 oz. (113 g)
NC98200-80	5 lbs. (2.3kg)

Medium-Soft Pink

NC98210-02	2 oz. (57 g)
NC98210-03	3 oz. (85 g)
NC98210-04	4 oz. (113 g)
NC98210-80	5 lbs. (2.3kg)

Medium Green

NC98220-02	2 oz. (57 g)
NC98220-03	3 oz. (85 g)
NC98220-04	4 oz. (113 g)
NC98220-80	5 lbs. (2.3kg)

Medium-Firm Blue

NC98230-02	2 oz. (57 g)
NC98230-03	3 oz. (85 g)
NC98230-04	4 oz. (113 g)
NC98230-80	5 lbs. (2.3kg)

Firm Purple

NC98240-02	2 oz. (57 g)
NC98240-03	3 oz. (85 g)
NC98240-04	4 oz. (113 g)
NC98240-80	5 lbs. (2.3kg)

Note: 3 oz. is our most popular size. 3 oz. (90 cc of volume) of Rainbow™ Exercise Putty is the ideal amount for large or small hands.

BEST VALUE

3 oz.cup

As low as

\$2.38

(60% OFF RETAIL)

Putty Containers

- Additional empty containers for exercise putties.
- Cups can hold up to 5 oz. (141 g) of putty.
- Sold in a pack of 10 cups.

NC51025 Putty Containers (Pack of 10)



To order, log in at www.ncmedical.com

MUSCLE AIDTAPE™

ARM YOUR MUSCLE **BEAT THE PAIN**

- Relieve muscle pain, soreness and cramping.
- Provide support for joints and muscles.
- Reduce swelling, inflammation and bruising.
- Increase blood circulation/flow.
- Water-resistant and breathable.
- Not made with natural rubber latex.



See product listings for available colors.



Available in:

Single rolls: 2" x 16' (5 cm x 5 m)

Bulk rolls: 2" x 105' (5 cm x 32 m)

Display Cartons	Color	Quantity
MT1642-MVBG	Beige	6 Rolls
MT1642-JDBK	Black	6 Rolls
MT1642-RPRD	Red	6 Rolls
MT1642-RPOR	Orange	6 Rolls
MT1642-MVPL	Purple	6 Rolls
MT1642-MVPK	Pink	6 Rolls
MT1642-JDNB	Navy Blue	6 Rolls
MT1642-JDBL	Blue	6 Rolls
MT1642-MVYW	Yellow	6 Rolls
MT1642-RPWH	White	6 Rolls
MT1642-RPUS	USA	6 Rolls

Bulk Rolls	Color
MT1052-BG-BULK	Beige
MT1052-BK-BULK	Black
MT1052-RD-BULK	Red
MT1052-NB-BULK	Navy Blue
MT1052-BL-BULK	Blue

As low as
\$6.57 /Roll
(40% OFF RETAIL)

**Strengthen Performance.
Relieve Muscle Pain.
Facilitate Rapid Recovery.**



Also Available in:
10" Long Precut Strips

2" x 10" (5 cm x 25.5 cm)

Precuts (20)

Color

MTPCBG-20	Beige
MTPCBK-20	Black
MTPCBL-20	Blue
MTPCUS-20	USA

Bulk Precuts (128)

Color

MTPCBG-128	Beige
MTPCBK-128	Black
MTPCUS-128	USA

**MUSCLE
AIDTAPE™**

Longer lasting, skin-friendly support.

To order, log in at www.ncmedical.com



Exercise Balls

Economical burst-resistant balls provide versatile exercise workouts.

- Tested to support up to 550 lbs. (250 kg).
- Helps to improve core stabilization by training the spinal column muscles as the body works to stay in balance.
- Balls are available individually in retail box with handle.
- Not intended for heavy lifting.
- Not made with natural rubber latex.

Ideal for clinic or home use.

Exercise Balls

	Color	Diameter
NC50100	Blue	45 cm (17 ¾")
NC50101	Red	55 cm (21 ½")
NC50102	Teal	65 cm (25 ½")
NC50103	Dark Blue	75 cm (29 ½")
NC50104	Silver	85 cm (33 ½")



Mini Ball

Strengthen key muscle groups, core health, endurance, posture and balance.


- High quality, 9" (23 cm) diameter therapy ball intended for exercising and strengthening the abdomen, back, hips, knees, shoulders, and neck.
- Anti-burst vinyl material deflates slowly if punctured.
- Designed to remain elastic and squishy when fully inflated.
- Supports up to 330 lbs. (150 kg).
- Not made with natural rubber latex.

As low as
\$5.37
 (40% OFF RETAIL)

NC50109 Mini Ball

Includes exercise poster.

To order, log in at www.ncmedical.com

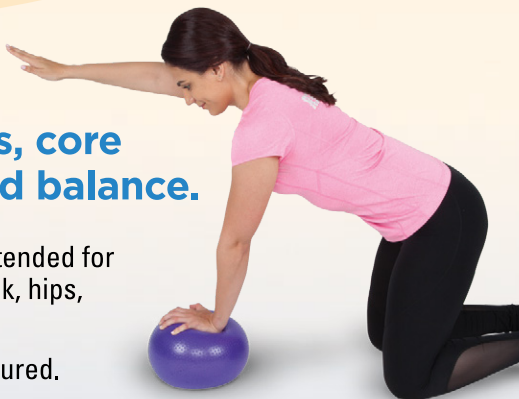

Stretch Strap

Safe and effective deep,
gradual stretching.



ISOTONICS
ISOMETRICS
PROLONGED STRETCH
INCREASED R.O.M.
FLEXIBILITY

- Designed with ten multi-positioned grip loops for use in deep, gradual stretching of the major muscle groups.
- Ideal for individuals seeking unassisted or assisted stretching exercises to relieve sore or achy muscles.
- Stretch more safely and effectively than is possible unaided.
- Grip loops can be used to set stretch improvement goals.
- Made of non-stretch webbing.
- Machine washable.
- Not made with natural rubber latex.

Retail ready packaging includes an instruction booklet for quick reference.

As low as
\$5.37
(40% OFF RETAIL)



Stretch Strap

	Model	Color	Length
NC84527	Standard	Purple	6' 4" (1.9 m) for people under 6' (1.8 m)
NC84528	Long	Blue	7' 5" (2.3 m) for people over 6' (1.8 m)

To order, log in at www.ncmedical.com

Common Trigger Points



BODYtool™

TRIGGER POINT SELF MASSAGER

**GREAT FOR: NECK, SHOULDERS,
BACK, ARMS, LEGS, TORSO & FEET!**

- Designed to target trigger points and muscle knots in the body to assist with temporary release of muscular tension.
- **Textured grip handles** allow better control for effectively targeting trigger points.
- Features **seven large knobs**, strategically placed to control pressure to the desired areas for trigger point release and breaking up of myofascial tissue.

NC69700 Body Tool™

Receive a
free retail display
when you buy
12 Body Tools™

As low as
\$23.97
(40% OFF RETAIL)



Easy-to-follow instruction
guide and poster included!



Retail
display rack
holds 12
Body Tool™
massagers.



**Body Tool™
Retail Stand
NC69700-RS**

To order, log in at **www.ncmedical.com**



Polyvinyl Cold Packs

Professional cold packs remain flexible
and comfortable down to **-4°F (-20°C)**.



NC70595

Oversize • 22 x 13 in. (56 x 33 cm)



NC70596

Neck Contour • 23 x 8 in. (58 x 20 cm)



NC70593

Half Size

9½ x 5½ in. (24 x 14 cm)



NC70594

Standard • 15 x 12 in. (38 x 30 cm)



- Designed for professional or at-home use.
- Effective relief of pain, discomfort or swelling after injury for up to 30 minutes.
- Conforms to body contours while providing soothing therapy.
- Filled with a premium, non-toxic gel that remains pliable and withstands freezing solid.
- Guaranteed not to split, crack or leak.
- One year manufacturer warranty.
- Not made with natural rubber latex.

As low as
\$5.96

(25% OFF RETAIL)

To order, log in at
www.ncmedical.com

Quality Products for Evaluation

Exacta™ Finger Pulse Oximeter

Accurate results in the clinic or at home.

- Dual color OLED screen displays SpO2, pulse rate waveform and pulse bar.
- Press the power button to easily change the direction of display.
- Anti-motion technology provides an accurate measurement.
- Resistant to interference by surrounding and man-made light.
- Measurement at low perfusion.
- Medical grade, FDA cleared.

NC87031 Finger Pulse Oximeter



Dimensions: 2.24" x 1.3" x 1.2"
(57 x 33 x 30 mm)

Weight: 1.9 oz. (54g)
with batteries.



Includes two "AAA" batteries with more than 18 hours operation.



Wrist strap included.

To order, log in at **www.ncmedical.com**

Quality Products for Evaluation

Exacta™ Devore® Pocket Goniometer

Accurate, economical and easy-to-read.

- Use to measure active and passive motion of the MP hand joints and the interphalangeal joints.
- Can be used on the volar surface of the fingers to measure hyperextension or laterally along the long axes of the joint.
- Measures from 0° to 180° in easy-to-read 1° increments.

NC70115 Devore® Pocket Goniometer



**EVERYDAY
EVERYONE™
PRICING**

**New
low price!**

\$8.95

To order, log in at
www.ncmedical.com

Wireless "On-the-Go" Technology

Intelligent iontophoresis.



Technologically advanced,
self-contained iontophoresis
drug delivery system.

- Fully self-contained.
- Intelligent biofeedback microprocessor.
- Hydrogel dispersive technology—no saline required.
- Patented "shut-off" technology.
- SmartPower™ LED.
- Detects and measures skin impedance to optimize drug delivery.



Active area: 11.4 cm²
Dosage: 40 or 80 mA*min
Skin Adhesive: Hypoallergenic
Maximum Buffer Capacity: Ag/AgCl Stabilizer
Fill Volume: 2.0cc
Smart Power: 6V with intelligent microprocessor
Packaging: 6 ActivaPatch® Kits, Alcohol Wipes,
Patient Information Cards
Treatment time: 2.5 Hours

**80 mA*min dose
in only 2.5 hours.**

**ActivaPatch® IntelliDose™ 2.5
NC89188**

* Study data on file

no saline needed.

To order, log in at www.ncmedical.com

ActivaPatch® IontoGo™ technology.

No external power source to connect, no leadwires,
no dose controller needed, power shut down feature and
no saline required.

ACTIVApatch®
IontoGo™ 4.0

Extended wear patch with an average drug delivery dose of 80 mA*min in approximately 4 hours.

4
Hours

ACTIVApatch®
IontoGo™ 12.0

Extended wear patch with an average drug delivery dose of 80 mA*min in approximately 12 hours.

12
Hours

**HOW SMART
IS YOUR
WIRELESS IONTO?**

Active area: 11.4 cm²
Dosage: 80 mA*min
Skin Adhesive: Hypoallergenic
Fill Volume: 2.0cc
Packaging: 6 IontoGo™ 4.0,
6 Alcohol Wipes
Treatment Time: 4 Hours

**NO
SALINE
REQUIRED**

**ActivaPatch®
IontoGo™ 4.0
NC89189**

**FULLY
SELF-
CONTAINED**

Active area: 11.4 cm²
Dosage: 80 mA*min
Skin Adhesive: Hypoallergenic
Fill Volume: 2.0cc
Packaging: 6 IontoGo™ 12.0,
6 Alcohol Wipes
Treatment Time: 12 Hours

**NO
SALINE
REQUIRED**

**ActivaPatch®
IontoGo™ 12.0
NC89191**

**FULLY
SELF-
CONTAINED**

Smarter Iontophoresis™



* Study data on file

To order, log in at **www.ncmedical.com**

Iontophoresis Delivery Kit



Variety of sizes and shapes offer great flexibility for iontophoresis treatments.

- Made of non-sensitizing, hypoallergenic, medical grade materials.
- Delivers ions with a pH balance up to 80 mA-min.
- Maximum strength, precise buffering with unique carbon- Ag/AgCl technology.
- The highly absorbent drug matrix has no messy hydrogel.
- Large ground electrode surface area lowers sensation and increases comfort.
- Each kit contains 12 individual treatment sessions that include: one activity electrode; one return electrode; one alcohol prep pad.



Advanced, technology-based manufacturing for a thinner, more pliable electrode.



Small
Square 1.5 cc
NC89250B



Medium
Square 2.5 cc
NC89251B



Large
Square 4.0 cc
NC89252B



Butterfly 2.0 cc
NC89255B

As low as
\$39.95
Any size

Not made with natural rubber latex.



High-contrast display indicates dose, time and current.

ActivaDose® II

Advantages

- More power than dual-channel devices for better penetration; up to 80 volts.
- Automatic current ramp-up and ramp-down for maximum comfort. Dose range from 0 to 80 mA-min with 0.1 mA-min increments.
- Time calculations are performed automatically and recalculated if dose and current settings are changed, or with pause/restart.

NC89186 **ActivaDose® II Controller**

To order, log in at **www.ncmedical.com**



Latex Free*

Eco-Stim™

Neurostimulation Electrodes



- Premium quality Carbon/Hydrogel for maximum comfort.
- Reusable up to 20 times when stored as recommended.
- Compatible with all NMS and TENS devices.
- Personalized packaging for patient data usage.
- FDA cleared for OTC purchase.
- 2mm DIN pigtail connector.
- Available in a variety of shapes and sizes in both foam or cloth

* Not made with natural rubber latex.

Cloth


Foam

Personalized packaging
for patient data usage.

As low as
\$1.38

(65% OFF RETAIL)

**4 electrodes
per pack**

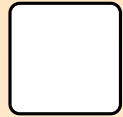


**2" x 2"
Square**

Foam: NC89292-1
Cloth: NC89293-1




**2" x 4"
Rectangle**

Foam: NC89292-3
Cloth: NC89293-3

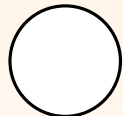

**2" diameter
Round**

Foam: NC89292-2
Cloth: NC89293-2

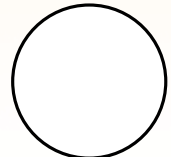

**2 3/4" diameter
Round**

Foam: NC89292-4
Cloth: NC89293-4




To order, log in at
www.ncmedical.com



redcoral®
PREMIUM NEEDLE SUPPLIERS

redcoral™ J-TYPE

ACUPUNCTURE NEEDLES

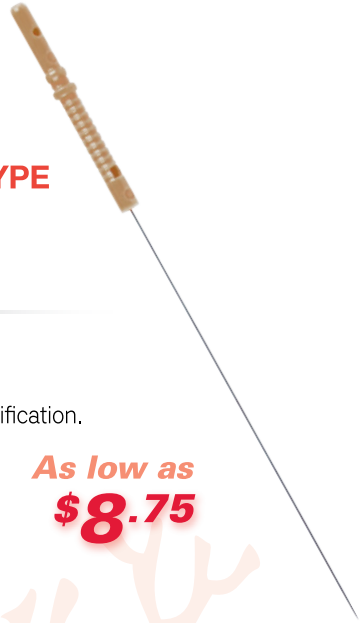

Premium quality needles for painless insertion.

- High grade German surgical steel.
- Parylene coated needle body for painless insertion.
- Color-coded molded plastic handles for easy size identification.
- No holding tab (time saving and no mess).
- One hand aim and insert application.
- 100 sterile needles per box.
- FDA certified.

As low as
\$8.75



Part #	Diameter	Length	Handle Color
RCJ1615	0.16 mm	15 mm	Red
RCJ1630	0.16 mm	30 mm	Red
RCJ1815	0.18 mm	15 mm	Ivory
RCJ1830	0.18 mm	30 mm	Ivory
RCJ1840	0.18 mm	40 mm	Ivory
RCJ2030	0.20 mm	30 mm	Blue
RCJ2040	0.20 mm	40 mm	Blue
RCJ2215	0.22 mm	15 mm	Pink
RCJ2230	0.22 mm	30 mm	Pink
RCJ2240	0.22 mm	40 mm	Pink
RCJ2530	0.25 mm	30 mm	Purple
RCJ2540	0.25 mm	40 mm	Purple
RCJ2550	0.25 mm	50 mm	Purple
RCJ3030	0.30 mm	30 mm	Brown
RCJ3040	0.30 mm	40 mm	Brown
RCJ3050	0.30 mm	50 mm	Brown
RCJ3060	0.30 mm	60 mm	Brown
RCJ3070	0.30 mm	70 mm	Brown
RCJ30100	0.30 mm	100 mm	Brown
RCJ3575	0.35 mm	75 mm	Brown
RCJ35100	0.35 mm	100 mm	Brown
RCJ40100	0.40 mm	100 mm	Brown

Stick with the best™

To order, log in at www.ncmedical.com


The Ideal Needle for Dry Needling **Just Got Better**

MYOTECH[®] DRY NEEDLES 2.0 ELITE

Produced in collaboration with healthcare practitioners and researchers specifically for dry needling.

- High grade German surgical steel.
- Parylene coated needle body for maximum comfort.
- Double clamped, 30% thicker aluminum handle improves dexterity and strength.
- Glue-free pink guide tube for efficient aim and tap insertions.*
- Microscopic grooves provide surface micro-channeling for near painless insertion.
- 100 sterile needles per box.
- FDA certified.

*75 mm-100 mm lengths feature an opaque guide tube.
135 mm length features clear guide tube with glue spot.

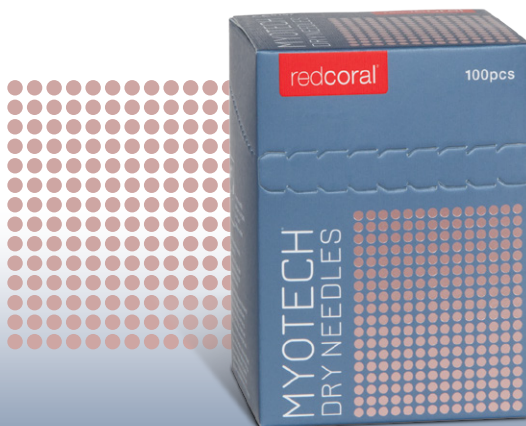

Part #	Diameter	Length
RCM2-1630	0.16 mm	30 mm
RCM2-1830	0.18 mm	30 mm
RCM2-1840	0.18 mm	40 mm
RCM2-2215	0.22 mm	15 mm
RCM2-2530	0.25 mm	30 mm
RCM2-2540	0.25 mm	40 mm
RCM2-2550	0.25 mm	50 mm
RCM2-2570	0.25 mm	70 mm
RCM2-3030	0.30 mm	30 mm
RCM2-3040	0.30 mm	40 mm
RCM2-3050	0.30 mm	50 mm
RCM2-3060	0.30 mm	60 mm
RCM2-3075	0.30 mm	75 mm
RCM2-30100	0.30 mm	100 mm
RCM2-40100	0.40 mm	100 mm
RCM2-45135	0.45 mm	135 mm

To order, log in at www.ncmedical.com



Better Tools, Better Outcomes™

Professional **PDN** Series & Specialist **SDN** Series

- Fully-lubricated needle body with medical grade, bio-accepting lubricant
- Surgical grade stainless steel gives shaft a firm feel at smaller diameter
- Conductive metal spring handle for ease of electrical stim application
- Shim needle placement (versus glue attachment)
- Precision tip designed for reduced insertion discomfort



PROFESSIONAL SERIES DRY NEEDLES

PDN Series

Sold in boxes of 100.



A premium quality and industry-specific product in the most common clinically used dry needle sizes for dry needling professionals.

	Diameter	Length
ID-PDN3030	.30 mm	30 mm
ID-PDN3040	.30 mm	40 mm
ID-PDN3050	.30 mm	50 mm
ID-PDN3060	.33 mm	60 mm

SPECIALIST SERIES DRY NEEDLES

SDN Series

Sold in boxes of 100.



Built for comfort, precision, and maneuverability when working in more sensitive body regions.

	Diameter	Length
ID-SDN2530	.25 mm	30 mm
ID-SDN2540	.25 mm	40 mm

To order, log in at www.ncmedical.com

CONVENIENT+PORTABLE+INFLATABLE



**EZ-UP
PILLOW™**

***Simple and comfortable solution
to elevate the legs or back.***

- Designed by a physical therapist as a portable, inexpensive and simple-to-use therapy pillow.
- Can be used to help relieve pressure, pain or swelling.
- Convenient and easy to use at home and while traveling.
- Clean easily with mild soap and water.
- Dimensions (inflated): 22" L x 22" W x 16" H (56 x 56 x 41 cm).
- Made of gray, non-phthalate PVC plastic.



*Space-saving
wedge deflates
for easy storage.*

Handle for easy carrying
Foot pump included

**EZ-UP Pillow™
NC80013**



To order, log in at **www.ncmedical.com**

CUSTOM ORTHOTICS

HEAT MOLDABLE & CUSTOMIZABLE, WITH 3 DENSITY OPTIONS TO CHOOSE FROM



Full
Length



3/4
Length

- Helps relieve pain in the heel, knee, and back.
- Combines durability with solid biomechanical support and control.
- Two lengths available: **Full** and **3/4**.

- **Low density** is ideal for patients that require special consideration such as seniors, and or patients with diabetes or arthritis.
- **Medium density** is ideal for patients with average weight and activity level and general orthotic needs.
- **High density** is ideal for athletes, high activity individuals, and bariatric patients.
- Sold in pairs.

Low Density Green

Part #	Size	Length	Men's Shoe Size	Women's Shoe Size
VAS71FG-XS	X-S	Full	3½ to 5	4½ to 6
VAS71RG-XS	X-S	3/4	3½ to 5	4½ to 6
VAS71FG-S	S	Full	5½ to 7	6½ to 8
VAS71RG-S	S	3/4	5½ to 7	6½ to 8
VAS71FG-M	M	Full	7½ to 9	8½ to 10
VAS71RG-M	M	3/4	7½ to 9	8½ to 10
VAS71FG-L	L	Full	9½ to 11	10½ to 12
VAS71RG-L	L	3/4	9½ to 11	10½ to 12
VAS71FG-XL	X-L	Full	11½ to 13	—
VAS71RG-XL	X-L	3/4	11½ to 13	—

Medium Density Blue

Part #	Size	Length	Men's Shoe Size	Women's Shoe Size
VAS71FB-XS	X-S	Full	3½ to 5	4½ to 6
VAS71RB-XS	X-S	3/4	3½ to 5	4½ to 6
VAS71FB-S	S	Full	5½ to 7	6½ to 8
VAS71RB-S	S	3/4	5½ to 7	6½ to 8
VAS71FB-M	M	Full	7½ to 9	8½ to 10
VAS71RB-M	M	3/4	7½ to 9	8½ to 10
VAS71FB-L	L	Full	9½ to 11	10½ to 12
VAS71RB-L	L	3/4	9½ to 11	10½ to 12
VAS71FB-XL	X-L	Full	11½ to 13	—
VAS71RB-XL	X-L	3/4	11½ to 13	—

High Density Red

Part #	Size	Men's Length	Women's Shoe Size	Shoe Size
VAS71FR-XS	X-S	Full	3½ to 5	4½ to 6
VAS71RR-XS	X-S	3/4	3½ to 5	4½ to 6
VAS71FR-S	S	Full	5½ to 7	6½ to 8
VAS71RR-S	S	3/4	5½ to 7	6½ to 8
VAS71FR-M	M	Full	7½ to 9	8½ to 10
VAS71RR-M	M	3/4	7½ to 9	8½ to 10
VAS71FR-L	L	Full	9½ to 11	10½ to 12
VAS71RR-L	L	3/4	9½ to 11	10½ to 12
VAS71FR-XL	X-L	Full	11½ to 13	—
VAS71RR-XL	X-L	3/4	11½ to 13	—
VAS71FR-XXL	XX-L	Full	13½ to 15	—

VASYLI + ARMSTRONG II

SENSITIVE FEET
SHEAR-REDUCING ORTHOTIC

IDEAL FOR
THE TREATMENT OF
SENSITIVE FEET
AND DIABETIC FEET



- Ideal for patients with plantar sensitivity, loss of plantar fat pad, arthritis, blistering, and peripheral neuropathy.
- Helps reduce shear forces and plantar ulcerations.
- Wider forefoot profile, making it suitable for extra wide footwear.
- Contains an antimicrobial top cover.
- Shock absorbing layer.
- Sold in pairs.

Part #	Size	Men's Shoe Size	Women's Shoe Size
VAS71VAII-XS	X-Small	3½ to 5	4½ to 6
VAS71VAII-S	Small	5½ to 7	6½ to 8
VAS71VAII-M	Medium	7½ to 9	8½ to 10
VAS71VAII-L	Large	9½ to 11	10½ to 12
VAS71VAII-XL	X-Large	11½ to 13	—
VAS71VAII-XXL	XX-Large	13½ to 15	—

VASYLI | Signature Series



VASYLI + DANANBERG

PROFESSIONAL ORTHOTIC DEVICE WITH
1ST RAY FUNCTION ENHANCEMENT FEATURE

IDEAL FOR THE
TREATMENT OF
LOWER BACK PAIN
AND PLANTAR
FASCIITIS



- Removable proximal and distal plugs that make treatment of 1st ray function easier.
- Tri-compound construction provides cushion and support.
- Orthotic is heat moldable, allowing for personalization of arch and rearfoot control.
- Helps relieve pain in the heel, knee, or back.
- Ideal for walking and everyday use.
- Sold in pairs.

Part #	Size	Men's Shoe Size	Women's Shoe Size
VAS71HDF-XS	X-Small	3½ to 5	4½ to 6
VAS71HDF-S	Small	5½ to 7	6½ to 8
VAS71HDF-M	Medium	7½ to 9	8½ to 10
VAS71HDF-L	Large	9½ to 11	10½ to 12
VAS71HDF-XL	X-Large	11½ to 13	—
VAS71HDF-XXL	XX-Large	13½ to 15	—

VASYLI | Signature Series



North Coast Medical, Inc.
780 Jarvis Drive, Suite 100
Morgan Hill, CA 95037

Our Company

For over 40 years, North Coast Medical has maintained an outstanding reputation in the global rehabilitation market. We are more than just a wholesale distributor and manufacturer of products, equipment & supplies. North Coast is a strategic supply chain partner, providing innovative opportunities for our customers and business partners. Now a second-generation family business, we continue to invest in our future while steadfastly adhering to strong ethics and integrity. We invite you to browse our website and discover what is more than expected at North Coast Medical.

Pricing Disclaimer


The "As low as" professional wholesale pricing referenced in this catalog is typically associated with a reasonable quantity. Wholesale pricing is advertised and offered to our healthcare customers. The "As low as" price is advertised with a reference to a discount from the Suggested Retail Price.

The products featured in this catalog are only a small sample of what you'll find on our website. Thousands of medical and rehabilitative products are available at www.ncmedical.com - log in to view your professional pricing. Orders can also be placed by phone at **800.821.9319**.

FREE SHIPPING

On e-com orders over \$49.00

*Available for
qualified accounts.*

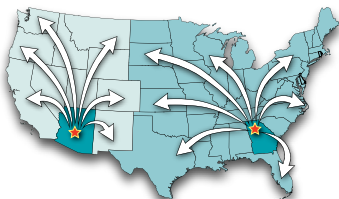

*Most orders delivered in
1-3 days anywhere in the U.S.*

 **North Coast**
MEDICAL & REHABILITATION PRODUCTS
More Than Expected™

NCC5550



7 68627 00559 3