
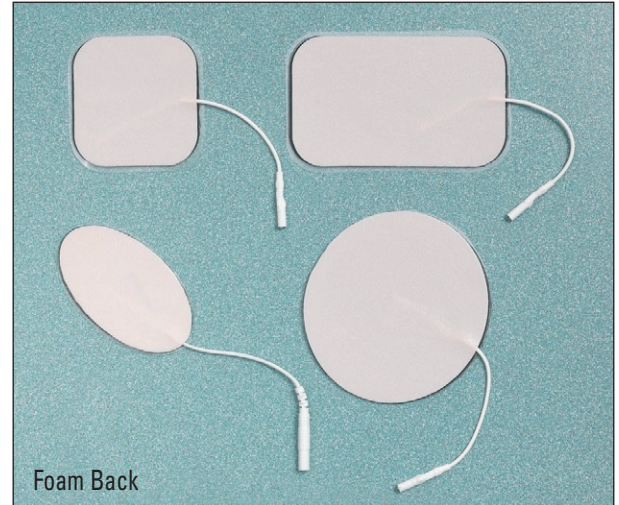



# Multi-use Electrodes

Comfort and economy in a variety of shapes and sizes

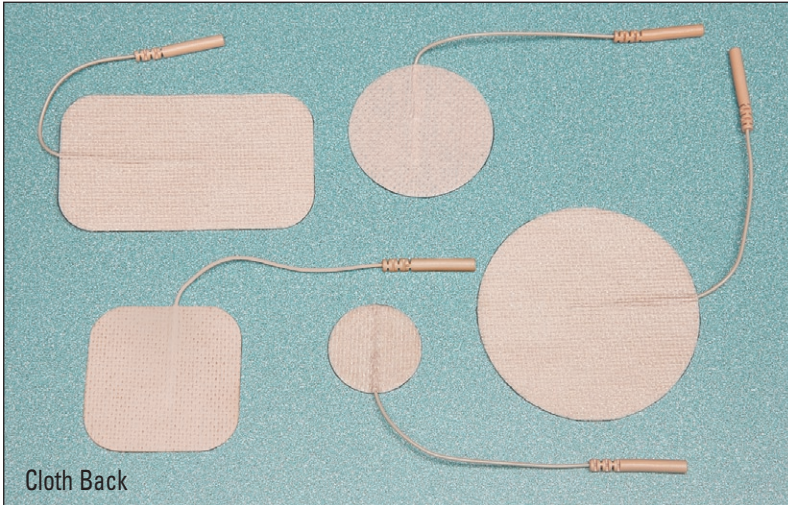

### North Coast™ Multi-Use Electrodes Provide Optimum Current Distribution

- Enhance patient comfort during stimulation.
- Disperse current evenly across the electrode.
- Eliminate stinging, edge biting and hot spots.

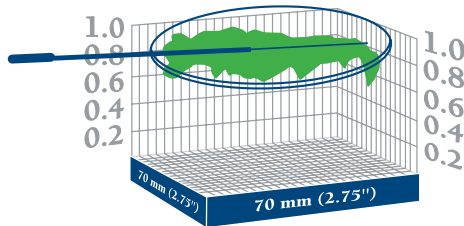

Published test methodology available upon request.

### MultiStick® patented two-layer adhesive gel

- Delivers uncompromised gel adhesion and performance.
- Eliminates performance problems associated with single-layer gels.
- Optimized for bonding to the electrode.
- Designed for multiple applications to the skin.



MultiStick® is a registered trademark of Axelgaard Manufacturing Co., Ltd.

## TENS and muscle stimulation electrodes with MultiStick® Gel.

These solid gel electrodes demonstrate exceptional conformity while maintaining superior conductive capabilities. Features dual-layer adhesive gel. Latex free. Each flexible pigtail contains 1,000 strands of stainless steel, which is placed precisely against the conductive material. This high-performance electrode is backed by cloth or foam.

### Foam Back

NC89281	Round	1¼"	(3.2cm)	Pkg. of 4
NC89282	Round	2"	(5.1cm)	Pkg. of 4
NC89283	Round	2¾"	(7.0cm)	Pkg. of 4
NC89284	Square	2" x 2"	(5.1 x 5.1cm)	Pkg. of 4
NC89285	Rectangle	2" x 3½"	(5.1 x 8.9cm)	Pkg. of 4
NC89287	Rectangle	3" x 5"	(7.6 x 13cm)	Pkg. of 2
NC89278	Oval	1½" x 2½"	(3.8 x 6.4cm)	Pkg. of 4
NC89279	Oval	2" x 4"	(5.1 x 10cm)	Pkg. of 4

### Cloth Back

NC89231	Round	1¼"	(3.2cm)	Pkg. of 4
NC89232	Round	2"	(5.1cm)	Pkg. of 4
NC89233	Round	2¾"	(7.0cm)	Pkg. of 4
NC89234	Square	2" x 2"	(5.1 x 5.1cm)	Pkg. of 4
NC89235	Rectangle	2" x 3½"	(5.1 x 8.9cm)	Pkg. of 4
NC89209	Rectangle	3" x 5"	(7.6 x 13cm)	Pkg. of 2



North Coast Medical, Inc.  
8100 Camino Arroyo • Gilroy, CA 95020  
Local/Int'l: 408-776-5000 • Toll-Free: 800-821-9319  
Toll-Free Fax: 877-213-9300 • www.ncmedical.com